

Power wheelchairs

# F5 Corpus VS




permobil

# Stand for more

The F5 Corpus VS that you know and trust is now available with Group 3 and Group 4 base options, making it more accessible than ever to pursue the functional and medical benefits of power standing.

	Group 3	Group 4
<b>Best-in-class base performance</b>	✓	✓
<b>Fully customizable standing sequences</b>	✓	✓
<b>Articulating power footplate</b>	✓	✓
<b>Functional drive speeds in standing</b>	✓	✓
<b>Full range of power tilt, recline, &amp; legrest elevation</b>	✓	✓
<b>14" of ActiveHeight™ seat elevation</b>	✓	✓
<b>Optional no charge accessories:</b> <i>LED light kit, USB charger, MyPermobil</i>	✓	✓
<b>Max speed of 7.5 mph</b>	✗	✓
<b>2-year warranty</b>	✗	✓
<b>Preventative maintenance plan</b>	✗	✓
<b>Hydrographic color options included</b>	✗	✓
<b>One-time shroud color change included</b>	✗	✓

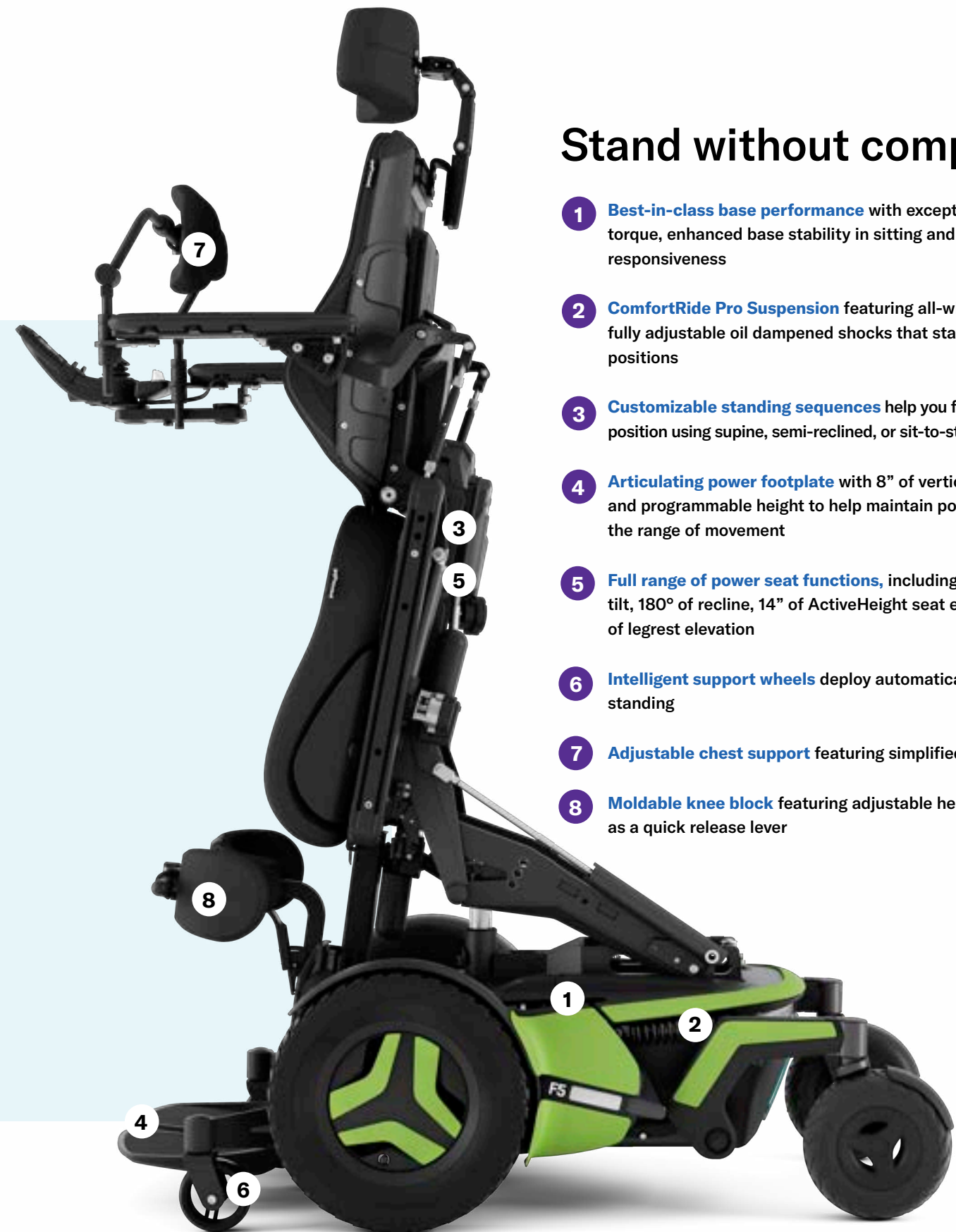

# Stand without compromise

- 1 Best-in-class base performance** with exceptional motor torque, enhanced base stability in sitting and standing & chair responsiveness
- 2 ComfortRide Pro Suspension** featuring all-wheel independent fully adjustable oil dampened shocks that stay active in all positions
- 3 Customizable standing sequences** help you find the perfect position using supine, semi-reclined, or sit-to-stand techniques
- 4 Articulating power footplate** with 8" of vertical foot plate travel and programmable height to help maintain positioning throughout the range of movement
- 5 Full range of power seat functions**, including: 50° of posterior tilt, 180° of recline, 14" of ActiveHeight seat elevation, and 180° of legrest elevation
- 6 Intelligent support wheels** deploy automatically for stability in standing
- 7 Adjustable chest support** featuring simplified adjustment points
- 8 Moldable knee block** featuring adjustable height & depth as well as a quick release lever

# Your outcome is our top priority

In addition to the superior power standing provided by the F5 Corpus VS, Permobil is proud to offer industry-leading clinical and funding resources. We want to ensure you and your healthcare team have what it takes to achieve an optimal outcome.



Scan the QR code to learn more



## Specifications

<b>Maximum user weight</b>	300 lb	<b>Drive wheel</b>	14"
<b>Maximum speed</b>	5.8 mph, Group 3 7.5 mph, Group 4	<b>Casters</b>	8"
<b>Range</b>	20 mi <sup>1</sup>	<b>Tilt options</b>	Posterior: 0° – 50° Standing
<b>Base width</b>	25.5"	<b>Legrest angle - power</b>	85° – 180°
<b>Base length</b>	43"	<b>Backrest angle - power recline</b>	85° – 180°
<b>Base turning radius</b>	30"	<b>Backrest height</b>	23" – 28" (1" incr.)
<b>Ground clearance</b>	3"	<b>Seat depth</b>	14" – 22" (1" incr.)
<b>Recommended battery type</b>	Group 24	<b>Seat width</b>	17" – 23" (2" incr.)
<b>Weight including batteries</b>	432 lb	<i>*Actual driving range will vary based on driving and battery conditions.</i>	
<b>Drive electronics</b>	R-Net 120A	<i><sup>1</sup>Measured without cushion.   Specifications may vary based upon configuration.</i>	
<b>Seat-to-floor height</b>	17.5", 18.5", 19.5" <sup>†</sup>	<i>Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and includes base, batteries, power seat functions and seating system without accessories.</i>	
<b>Seat elevation</b>	14"		